

FLA OFFERS THE MOST UNIQUE LINE OF NIGHT SPLINT AFO PRODUCTS AVAILABLE:

- Indications:
- Arthritis
 - Achilles Tendinitis
 - Plantar Fasciitis
 - Calf Muscle Cramps
 - Equinus (Muscle Tightening)
 - Heel & Arch Pain
 - Achilles Paratenonitis
 - Post Surgical Treatment of the Foot

HealWell® AFO Night Wrap

MODEL: 58-600

For New or Mild Conditions

- Extremely comfortable low profile with countervailing tension
- Ideal when neutral or slight dorsiflexion is desired
- Very comfortable for elderly arthritic patients

[More information on page 14](#)



- Indications:
- Plantar Fasciitis
 - Achilles Tendinitis
 - Equinus (Muscle Tightening)
 - Foot Drop
 - Calcaneal Apophysitis
 - Post Surgical Treatment of the Foot and Ankle

**HealWell® CUB™
Plantar Fasciitis Night Splint
Resting Comfort Slipper™**

MODEL: 58-500/501

For Moderate to Severe Conditions

- Especially useful when a controlled setting of dorsiflexion is needed or patient compliance is essential
- Half the bulk and twice the comfort

[More information on page 15](#)



- Indications:
- Plantar Fasciitis
 - Achilles Tendinitis
 - Equinus (Muscle Tightening)
 - Foot Drop
 - Non-Fixed Contractures
 - Pronation Syndromes

HealWell® Plantar Fasciitis Night Splint

MODEL: 58-409

For Moderate to Severe Conditions

- Especially useful for larger patients, those needing greater dorsiflexion or more rigid support, or patients with plantar flexion power

[More information on page 16](#)



FLA has compiled a monograph of abstracts from leading publications that demonstrate the effectiveness of night splints in the treatment of Plantar Fasciitis.

....Proven Plantar Fasciitis Relief!


Nearly 70% of the 2.5 million new cases of Plantar Fasciitis each year can be easily resolved...

Study after study demonstrates the effectiveness of night splints. Low cost and easy, their use for just 5 or 6 weeks can, in most cases, eliminate symptoms and speed healing.

A Conservative, Economical First Step...

...Saves you time and effort for those difficult cases!

To request a copy of this monograph, please contact Customer Service.

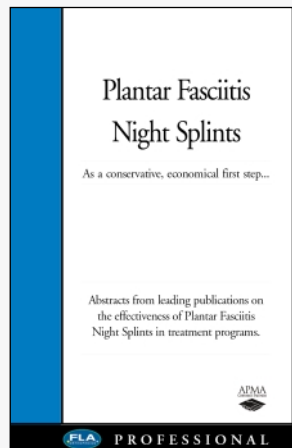


CALL FOR A COPY!



AFO Night Wrap

MODEL: 58-600

PATENTS
PENDING

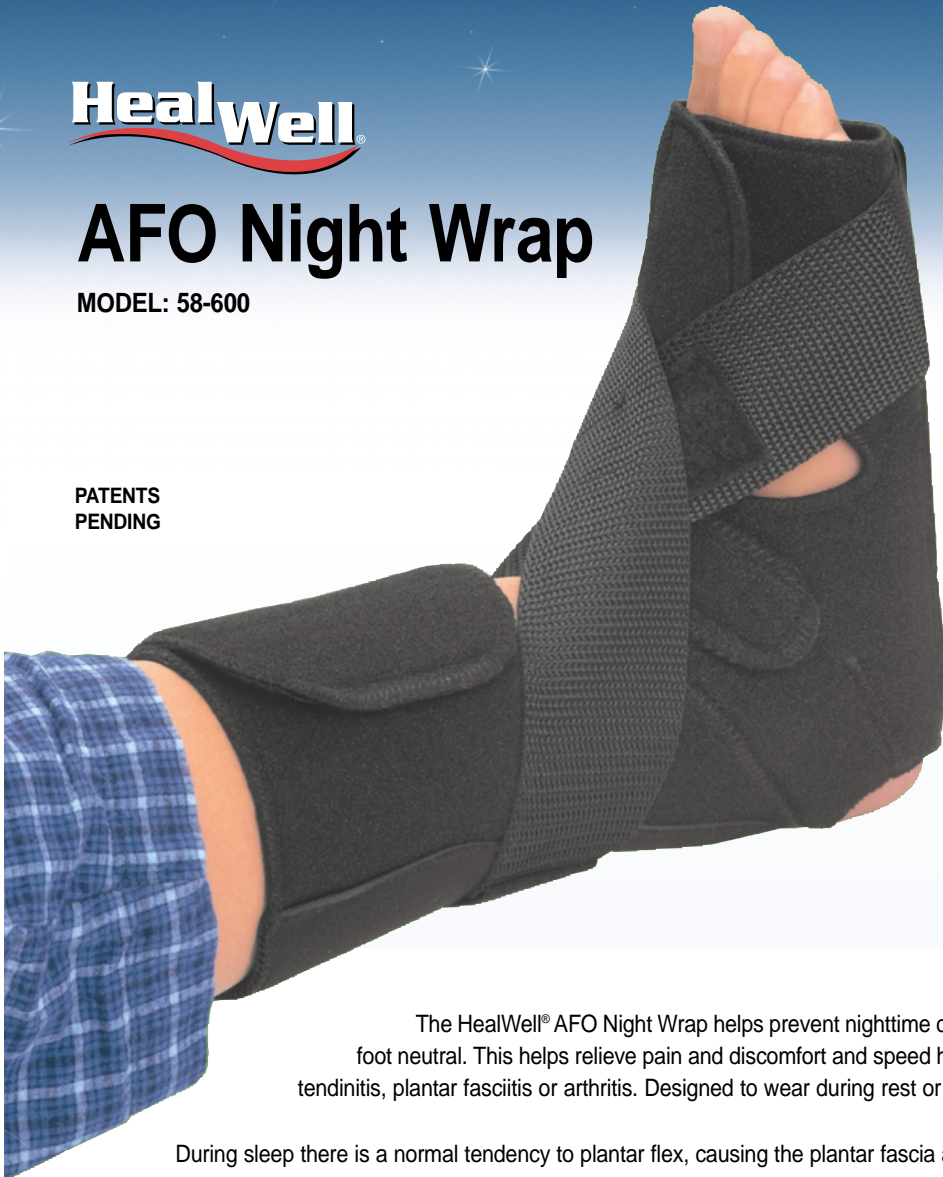

For New or Mild Conditions

- Extremely comfortable low profile with countervailing tension
- Ideal when neutral or slight dorsiflexion is desired
- Very comfortable for elderly arthritic patients

Indications:

- Arthritis
- Achilles Tendinitis
- Plantar Fasciitis
- Calf Muscle Cramps
- Equinus (Muscle Tightening)
- Heel & Arch Pain
- Achilles Paratenonitis
- Post Surgical Treatment of the Foot

The HealWell® AFO Night Wrap helps prevent nighttime contracture with countervailing tension to hold the foot neutral. This helps relieve pain and discomfort and speed healing in new or mild conditions such as Achilles tendinitis, plantar fasciitis or arthritis. Designed to wear during rest or sleep.

During sleep there is a normal tendency to plantar flex, causing the plantar fascia and Achilles tendon to contract from their normal active length. Nighttime contracture can result in significant heel and foot pain and discomfort for those with Achilles tendinitis, plantar fasciitis or arthritis. This simple device is designed to control nighttime contracture to treat and relieve discomfort in new or mild cases of these conditions at their onset to help prevent further injury. Color: Black.

Features:

- **SOFT LAMINATED FOAM** - Breathable, lightweight material is comfortable for night-time wear with a smooth Lycra lining.
- **SLIM, LOW PROFILE** - No bulk, compact wrap design is easy to sleep in.
- **NON-SLIP PAD** - For safe, occasional walking.
- **COMPRESSION FORMED BACK** - Holds the wrap in place and conforms around the ankle and heel for improved foot positioning.
- **DUAL TENSION STRAPS** - Countervailing tension straps help hold the foot in position and prevent plantar flexion. Straps criss-cross for less bulk.
- **PLANTAR SURFACE SUPPORT** - Tension straps work with a low profile, rigid stay on the plantar surface, functioning as a lever to maintain foot position.

**IDEAL FOR
ELDERLY ARTHRITIC
PATIENTS**

Fits Right or Left

ORDER #	SIZE	FITS	
		WOMEN'S SHOE	MEN'S SHOE
58-6004	Small/Medium	6 - 9	up to 8
58-6006	Large/X-Large	9½ and up	8½ to 11½

For shoe size greater than 11½ it is recommended to use the HealWell® Cub™ or Standard Night Splint