

SAFE-T-SPORT NEOPRENE SHOULDER SUPPORT

Model No: 16-500

Sizes: Small -Large

Color: Black

SIZING:

Measure chest circumference.

Size	Fits
Small	Up to 40"
Medium	41 - 50"
Large	51 - 60"

CAUTION: This product contains neoprene. Contact of neoprene material with skin may cause irritation. Persons with a known susceptibility to dermatitis or with broken skin should avoid contact with this material. If skin irritation develops from contact with this material, discontinue use immediately and consult a physician. Should not be worn by those with phlebitis or circulatory disorders. Not recommended for use during sleep.

Consult a physician if pain, swelling or unusual reactions are experienced. This device is only one element in the overall treatment program. There is no guarantee that injury will be cured or prevented through use of this product. Device should be used under the supervision of a medical professional.

The Safe-T-Sport® Neoprene Shoulder Support stabilizes, supports, and protects the shoulder to prevent reinjury, while still allowing for flexibility. The sports neoprene material gives comfortable compression and four-way stretch for a proper fit and comfort. Neoprene retains the body's natural heat for therapeutic and soothing warmth to the shoulder joints and surrounding muscles. Adjustable chest and arm closures allow for a custom fit. Built-in pocket holds a reusable hot/cold gel pack (sold separately) in place over the tender joint. Easy to wear under clothing. For correct sizing, measure around the chest. Color: Black.






- Protects and supports the shoulder while allowing for flexibility
- Sports neoprene with a soft lining surrounds the shoulder for therapeutic warmth and comfortable compression
- Built-in pocket holds a reusable hot/cold gel pack in place over the tender joint (gel packs sold separately)
- Use an optional thermal gel pack to provide hot or cold therapy
- Adjustable chest and arm closures
- Inconspicuous wear - easy to wear under clothing

Wear Instructions: This shoulder support should be worn against the skin under loose fitting clothing. Persons with acute injuries may require assistance with application and removal.

1. Unfasten the closure on the chest portion of the support (opposite of arm closure).
2. Insert arm of affected shoulder into the arm sleeve of the support. Pull support up over affected shoulder.
3. Secure chest closure under the opposite shoulder.
4. Readjust arm closure and both chest closures on either side to get the best fit if necessary.
5. Product when properly applied.

Hot/Cold Pack Use: The use of a hot/cold gel pack is optional. Follow instructions on the gel pack for hot or cold use and insert into pocket over shoulder. Recommended gel pack size is 6x10". Gel packs sold separately.

Care Instructions: Fasten all closures. Hand wash in cold water with mild soap. Wrap in a towel to remove excess water (do not twist or wring). Air dry away from heat. Do not tumble dry. Do not bleach. Wash frequently if worn during activities, such as work and sports.

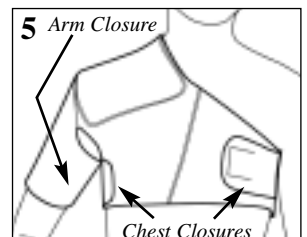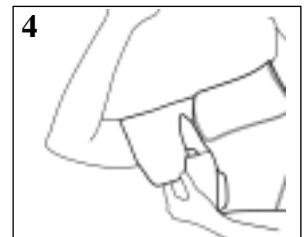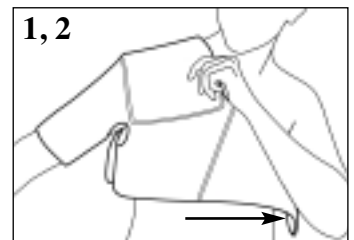

800-327-4110



Website: www.flaorthopedics.com

©2007 FLA ORTHOPEDICS, INC.