

POSTURA[®]

GELCELL™ CONTOUR CUSHION

*For Positioning,
Skin Protection
& Pressure Relief*

Models:

C3410	16 x 16 x 4"
C3450	18 x 16 x 4"
C3480	20 x 18 x 4"

Color: Black

FEATURES:

- **Anatomically Positioned Gel Insert** for maximum pressure relief under ischial tuberosities and coccyx
- **Dual Firmness Foam/Gel Design** with low durometer viscoelastic gel insert for maximum pressure relief and skin protection
- **Positioning Contours** – Durable, high resilience (HR) foam base is formed with leg troughs and thigh supports for optimal leg positioning. Provides a comfortable sitting position with lower extremity control and stability.
- **Fluid Proof Stretch Cover** – Removable, washable and wipeable
- **Cover is Breathable** for long lasting comfort. The loose fit allows the pelvis to “sink” farther into the cushion rather than restrict conformity to the user.
- **Cover has Soft, Shear Reducing Surface Fabric**
- **Cover meets CAL 117 Flammability Standard**
- **Handle** for easy carrying
- **Non-Slip Bottom** holds cushion securely in place to prevent accidental slippage or shifting



INDICATIONS:

For extra seated comfort and support. For individuals at medium to high risk of developing pressure ulcers or who need or prefer slight leg positioning for proper posture and added comfort. For those individuals who use a wheelchair for mobility but can ambulate and shift weight independently.

CARE INSTRUCTIONS:

Maintain by simply wiping down the fluid proof cover with warm water and mild soap. Air or towel-dry only. Machine wash removable cover **ONLY** in cold water with mild soap on gentle cycle. Air dry away from heat.

CUSHION DIRECTION:

Place the cushion in the desired seat with the Medial Thigh Support to the front of the chair and the gel insert at the back. With the cushion in proper place, sit in usual seated position with arm and foot rests properly placed.

SIZES:


Model	Width x Depth x Thickness
C3410	16" x 16" x 4"
C3450	18" x 16" x 4"
C3480	20" x 18" x 4"

Replacement covers available.

Cover Material:

Top: 90% Nylon, 10% Spandex
with Polyurethane Fluid Proof Film

Bottom : Coated Polyester

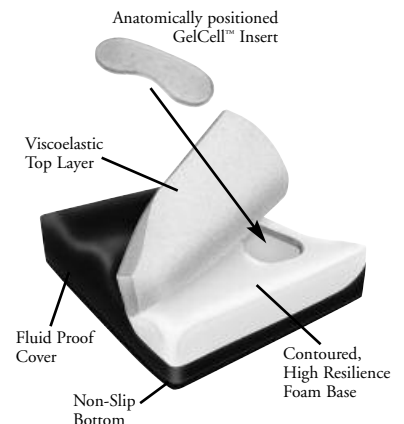

800-327-4110



Website: www.flaorthopedics.com

©2004 FLA ORTHOPEDICS, INC.