

# POSTURA®

## GELFOAM™ CUSHION

*For Skin Protection  
& Pressure Relief*

### Models:

*Viscoelastic Foam Top*

**C2335** 16 x 16 x 3.5"

**C2355** 18 x 16 x 3.5"

**C2375** 20 x 18 x 3.5"

*Convuluted Foam Top*

**C2457** 18 x 16 x 4"

**Color:** Black

### FEATURES:

- **Dual Firmness Foam/Gel Combination** for improved support and comfort
- **Pressure Relieving Low Durometer Viscoelastic Gel** in a continuous layer for superior skin protection and improved pressure/weight distribution.
- **Continuous Viscoelastic Gel Layer** will not displace or “bottom-out” like fluid gel bladders
- **Fluid Proof Cover** – Removable, washable and wipeable
- **Cover has Soft, Shear Reducing Surface Fabric.** The loose fit allows the pelvis to “sink” farther into the cushion rather than restrict conformity to the user.
- **Cover meets CAL 117 Flammability Standard** for greater protection
- **Handle** for easy carrying
- **Non-Slip Bottom** holds cushion securely in place to prevent accidental slippage or shifting



### Different Variations Available:

MODELS C2335, C2355, C2375: **Tri-Layer with Viscoelastic Foam and Viscoelastic Gel Layer over a firm foam base.** Unlike inexpensive water-based gel packs, this quality viscoelastic gel will never leak or bottom out. Provides additional pressure relief, comfort and improved weight distribution. Slow recovery viscoelastic foam top layer.

MODEL C2457: **Tri-Layer with Convuluted Foam and Viscoelastic Gel Layer over a firm foam base.** Viscoelastic gel layer provides additional pressure relief, comfort and weight distribution. Convuluted foam helps reduce risk of pressure ulcers by reducing shear and friction while increasing weight distribution and air circulation.

### INDICATIONS:

For extra seated comfort and support. For individuals at low to medium risk of developing pressure ulcers. For those individuals who use a wheelchair for mobility but can ambulate and shift weight independently.

### CARE INSTRUCTIONS:

Spot wipe clean removable cover with damp cloth. Machine wash outer cover ONLY in cold water with mild soap on gentle cycle. Air dry away from heat.

### SIZES:

Model	Width x Depth x Thickness	Top Layer
C2335	16" x 16" x 3.5"	Viscoelastic Foam
C2355	18" x 16" x 3.5"	Viscoelastic Foam
C2375	20" x 16" x 3.5"	Viscoelastic Foam
C2457	18" x 16" x 4"	Convuluted Foam


Replacement covers available.

Cover Material:


Top: 100% Treated Polyester  
with Polyurethane Fluid Proof Film

Bottom: Coated Polyester

*Model C2335, 2355 & 2375  
Viscoelastic Foam Top*

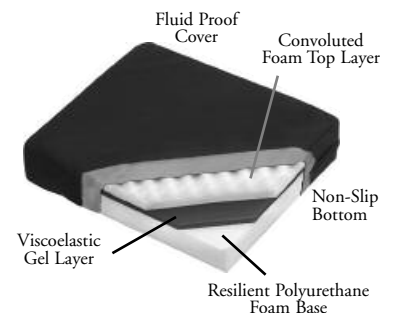

*Model C457  
Convuluted Foam Top*

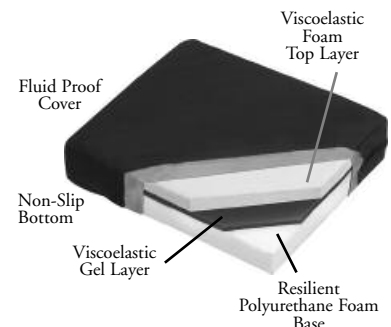

800-327-4110



Website: [www.flaorthopedics.com](http://www.flaorthopedics.com)

©2004 FLA ORTHOPEDICS, INC.